



## Putting On the L300 Go Functional Stimulation Cuff

- 1 While seated, slightly extend your lower leg.
- 2 Tilt the top of the L300 Go Cuff toward you, and then slide it up your leg until the locator rests below the kneecap. Lower the bottom of the L300 Go Cuff onto your leg.
- 3 Hold the L300 Go Cuff in place with your thumb or other hand on the EPG cradle. Pull the strap handle around the cradle.
- 4 Figure 4 shows the L300 Go Cuff correctly placed on the right leg.

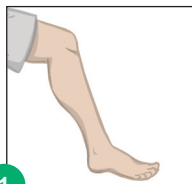

1

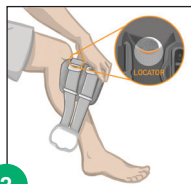

2










3

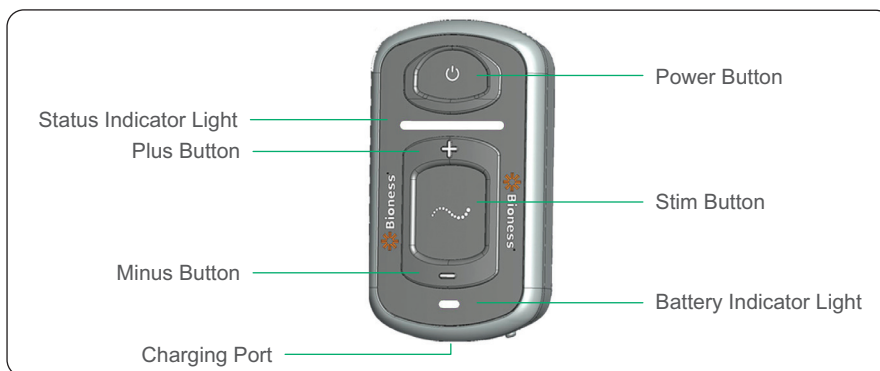

4

## Operating the L300 Go System







- 1 Press the  on/off button to turn on the EPG.
- 2 Select either gait or training mode:
  - To select gait mode, press the  mode button *briefly*.
  - To select training mode, press and *hold* the  mode button for 3 seconds. The EPG will acknowledge the new mode with audio and tactile feedback. To exit training mode, press and hold the  mode button for 3 seconds.
- 3 If necessary, adjust the stimulation intensity level by pressing the  plus or  minus button.



Charging the Batteries






Connect the system charger set to the EPG, then plug into a wall socket.



 (flashing)	EPG Status Indicator Bar Flashes Red along with Battery Indicatory Light Flashing Yellow	EPG Low Battery
 (flashing)	Flashes Yellow	Low Battery
 (flashing)	Flashes Green	Charging
 (flashing)	Solid Green	Fully Charged

EPG Displays










 (flashing)	Flashes Green	System is On
 (flashing)	Flashes Green	Stimulation is Off
 (flashing)	Flashes Yellow	Stimulation is On
 (flashing)	Flashes Red with Audible Alarm	Active Error
 (flashing)	Solid Red	EPG Malfunction







Control Unit Displays\*



\*Optional Component

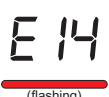	Solid	System is On
	Stim Wave Visible	System is in Gait/Training, Stimulation is enabled but Off
	Stim Wave Flashes	System is in Gait/Training, Stimulation is enabled and On
	Displays 0–9	Intensity Level
	Displays Gait Icon	Gait Mode is Enabled
	Displays Training Icon	Training Mode is Enabled
	Displays Audio Icon	Audio is On

System Error Displays

 (flashing)	Faulty Electrode Contact. Ensure electrodes are snapped in, EPG fully in cradle. Clear error alert by pressing  or  on EPG.		
 (flashing)	Low Battery Indicator on EPG. Charge system.		EPG Malfunction. Press and hold  button for 30 seconds then power on or call Technical Support.

For a complete list of all contraindications, warnings, adverse reactions, precautions, and care instructions please review the User's Guide.